



MULTITECH

021-22321127

09120783350



Combination Treatments For Body Slimming & Skin care

Skin Rejuvenation

- 1- Fractional RF Micro-Needling
- 2- Fractional RF No-Needling
- 3- RF Tightening Lift (bipolar)



Skin Care

- 1- Ulterasonic Face
- 2- Ulterasonic Eye
- 3- Iontophoresis
- 4- Microdermabration

Body Slimming

- 1- RF Lipolysis
- 2- Cavitation
- 3- Roller Vacuum

Fractional RF Micro-Needling




Microdermabrasion

Fractional RF No-Needling



RF Tightening Lift (bipolar)

Skin rejuvenation includes various treatments that aim to restore your skin from any damage. Skin damage could be a result of sun exposure, an underlying health condition, or a normal sign of aging.

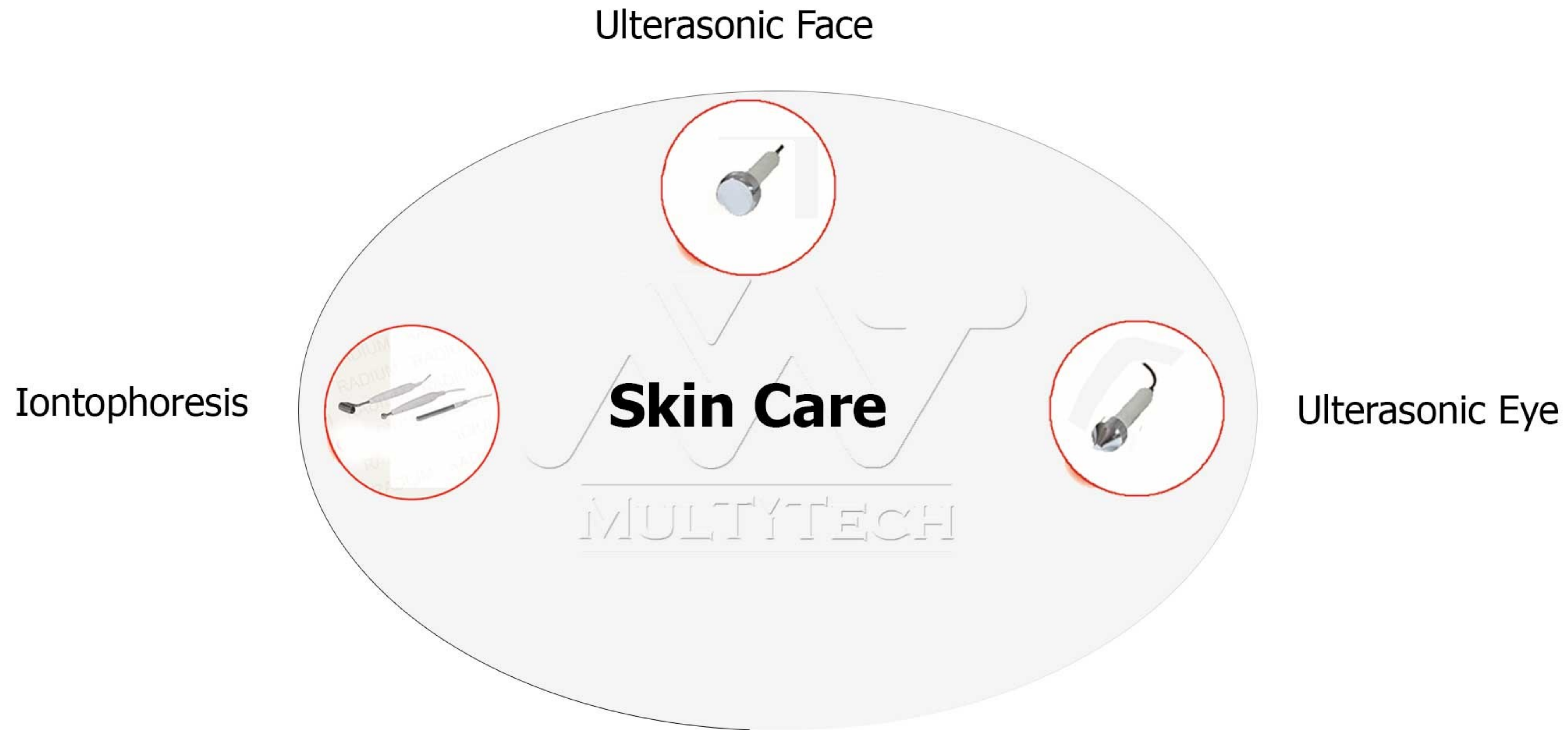

Skin care is the range of practices that support skin integrity, enhance its appearance and relieve skin conditions. They can include nutrition, avoidance of excessive sun exposure and appropriate use of emollients.



The term body slimming can also be used to refer to stimulation current treatment or cosmetic treatment. During stimulation current treatment, selective stimulation of nerves and muscles can be achieved locally

Fractional RF Micro-needling & No-needling

- *- Stimulation of fibroblasts and collagen production in the skin
- *- Reduction and removal of lines around the eyes and forehead
- *- Fixing drooping eyelids and face
- *- Aid in the absorption of topical skin medications
- *- Reducing open skin pores and helping to treat acne
- *- Reduction and removal of scars and stretch marks
- *- Neck Lift



RF Tightening Lift (bipolar)

- *-Fix sagging and laxity of the face
- *-Stimulating fibroblasts and rebuilding and building collagens
- *-Reduce and eliminate facial wrinkles
- *-Reduce and eliminate snoring

MULTITECH



Microdermabrasion

- *-Deep peeling of the skin
- *- Shrinking skin pores
- *- Reduction of pigment spots such as pregnancy spots and sunburn
- *- Scar removal


MULTITECH



Ulterasonic Face & Eye

- *- Increase blood circulation
- *- Help increase lymphatic flow
- *- Reducing puffiness and bruising around the eyes
- *- Reduce spots and freckles

MULTITECH

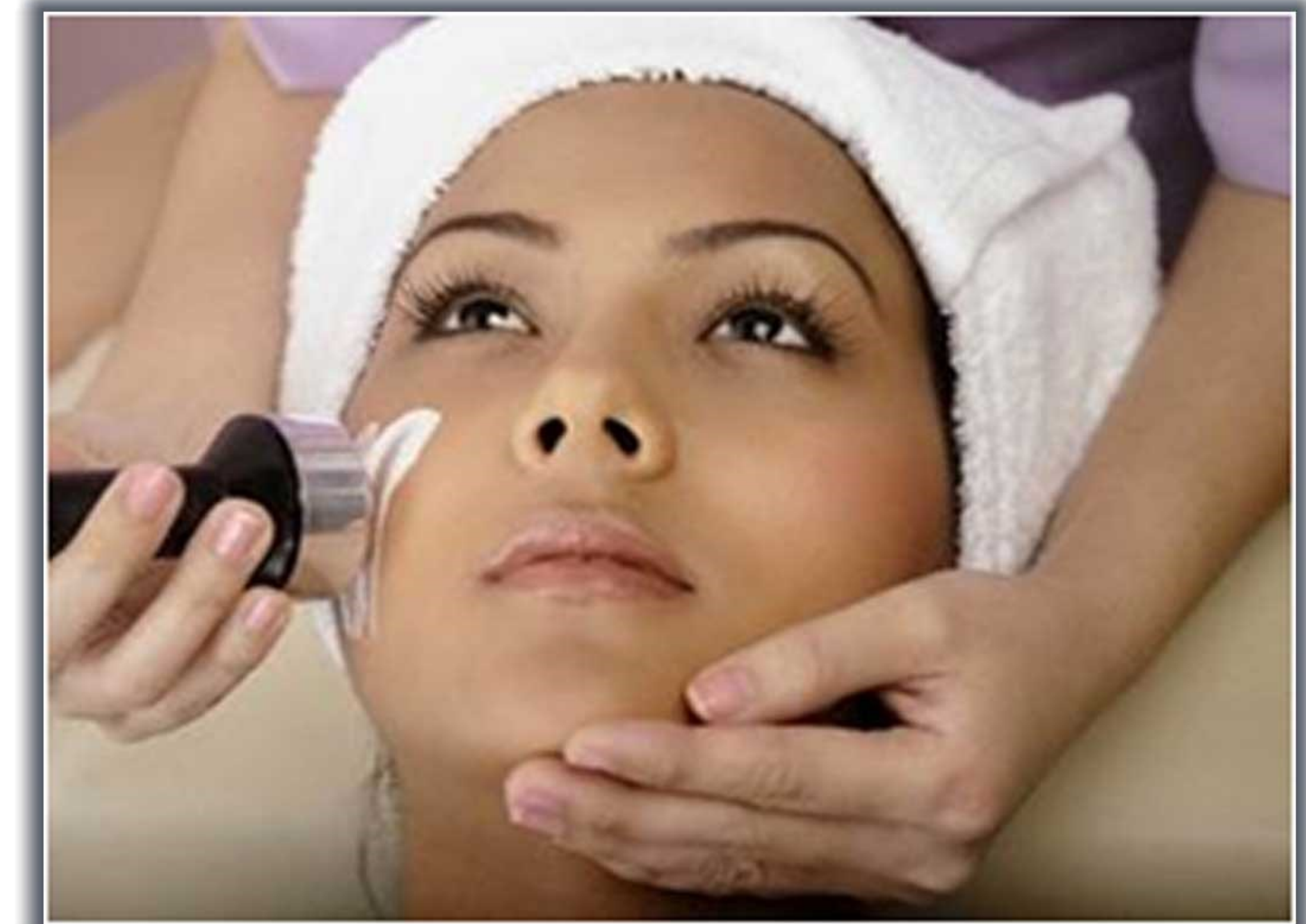

Iontophoresis

- *- Opening pores and deep cleansing of the skin
- *- Breakdown of hydrolipidic acid (fat)
- *- Stimulating the secretion of sebaceous glands and vessels
- *- Firmness and improvement of skin texture
- *- Hydroderma or skin hydration
- *- Adsorption of effective ionizable substances
(anti-stain, anti-wrinkle)



Cavitation

- *-Effective treatment of cellulite
- *-Removal of local fat and size reduction
- *- Fitness

MULTITECH



RF Lipolysis

- *-Removing local fat and reducing the size of the abdomen and sides
- *- Body skin lift
- *- Fitness
- *- Fix the squeamishness

MULTITECH



Roller Vacuum

- *- Increased capillary blood flow
- *- Removal of waste water and toxic substances from the site by lymphatic drainage
- *- Increase metabolism through oxygenation
- *- Reducing the size of fat cells
- *- Helps in the effective treatment of cellulite

MULTITECH



Technical Specifications

Combination Products System

Class IIa , type BF for (Mode RF tightening Lifting – Vacuum), type B for (Microdermabrasion)

Unit Type: Beauty Skin

Input Voltage: 220Vac

Mains Frequency: 50/60Hz

Max. Power Output: 25 Watt at 100 ohms (RF tightening Lifting)

Max. Power Output: 0.18Watt at 50 ohms (RF tightening Lifting)

Current Consumption with Max. R.F. Power: 0.5A

Output Frequency RF: 1MHz

Output Vacuum Pressure Microdermabrasion: 0.6Bar

Mass Of the Device: 30Kg

After and Before





After and Before




Microdermabrasion

Microdermabrasion is a cosmetic procedure that uses fine crystals and a vacuum to remove dead skin cells. It can be used on the face, neck, chest, back and hands.

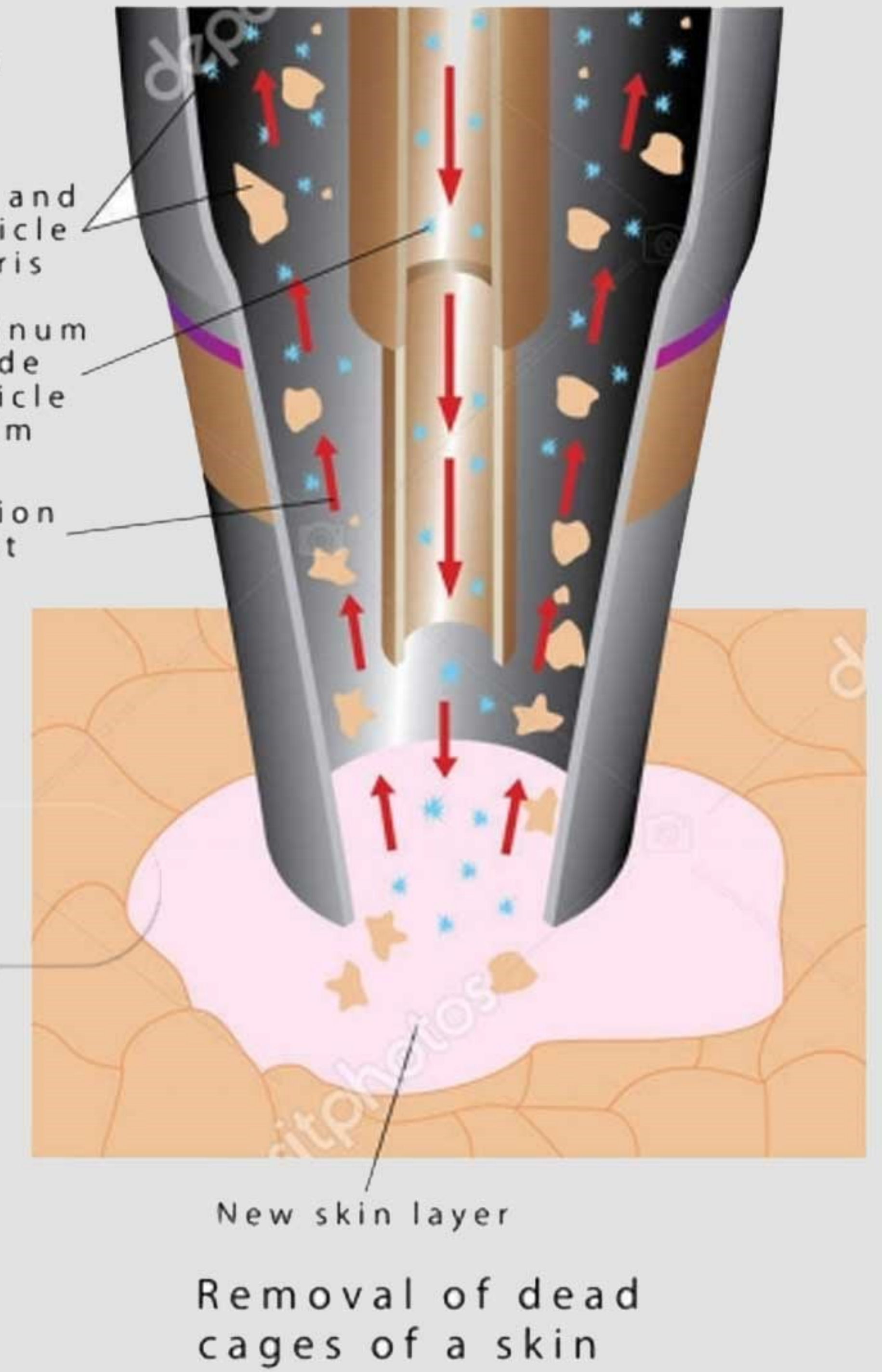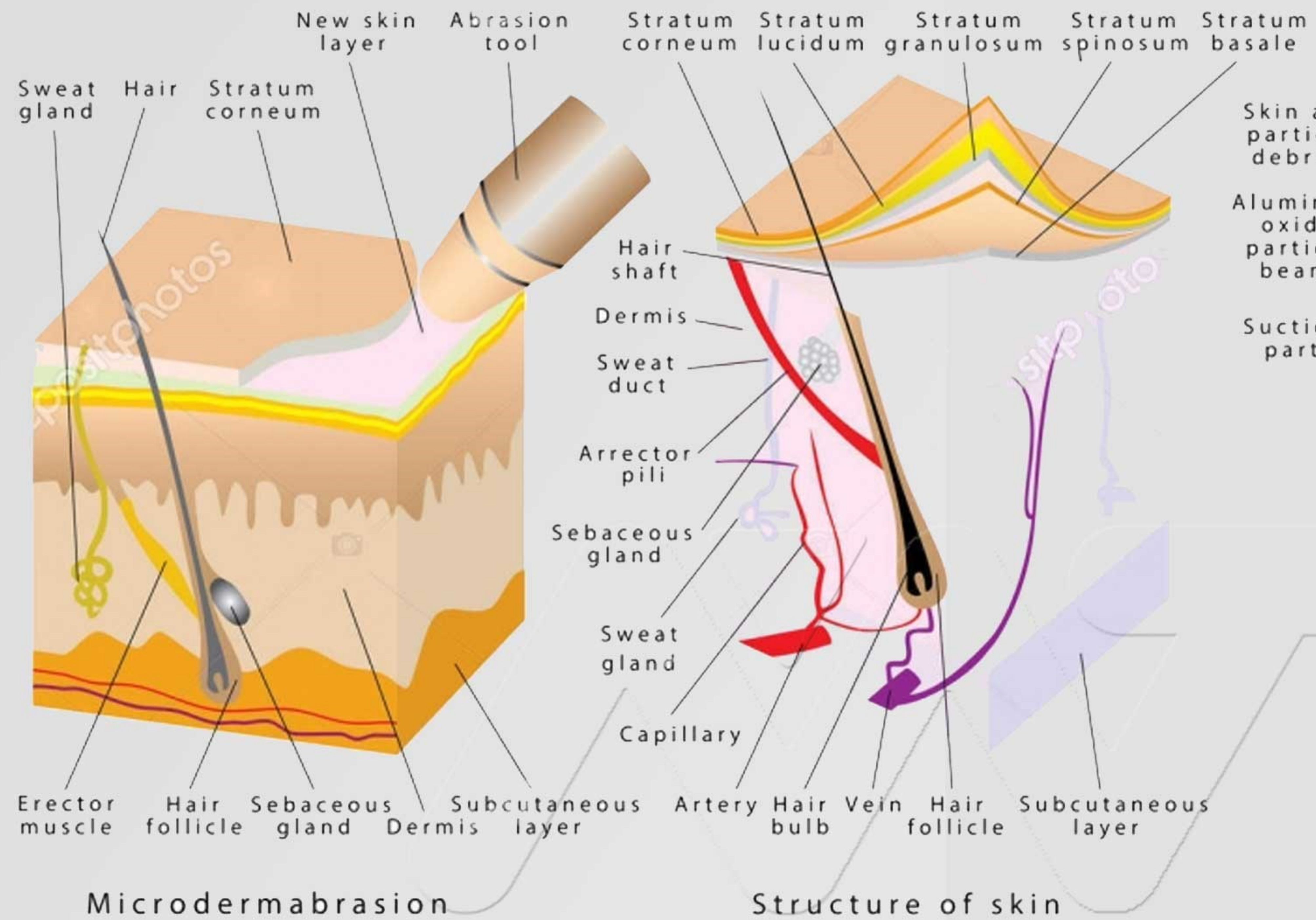


Fractional RF Micro-Needling & No-Needling


Fractional RF microneedling is a minimally invasive technology that involves the use of an array of very fine sterilized needles to create microscopic wounds at various predetermined depths in the dermis, without affecting the epidermis.

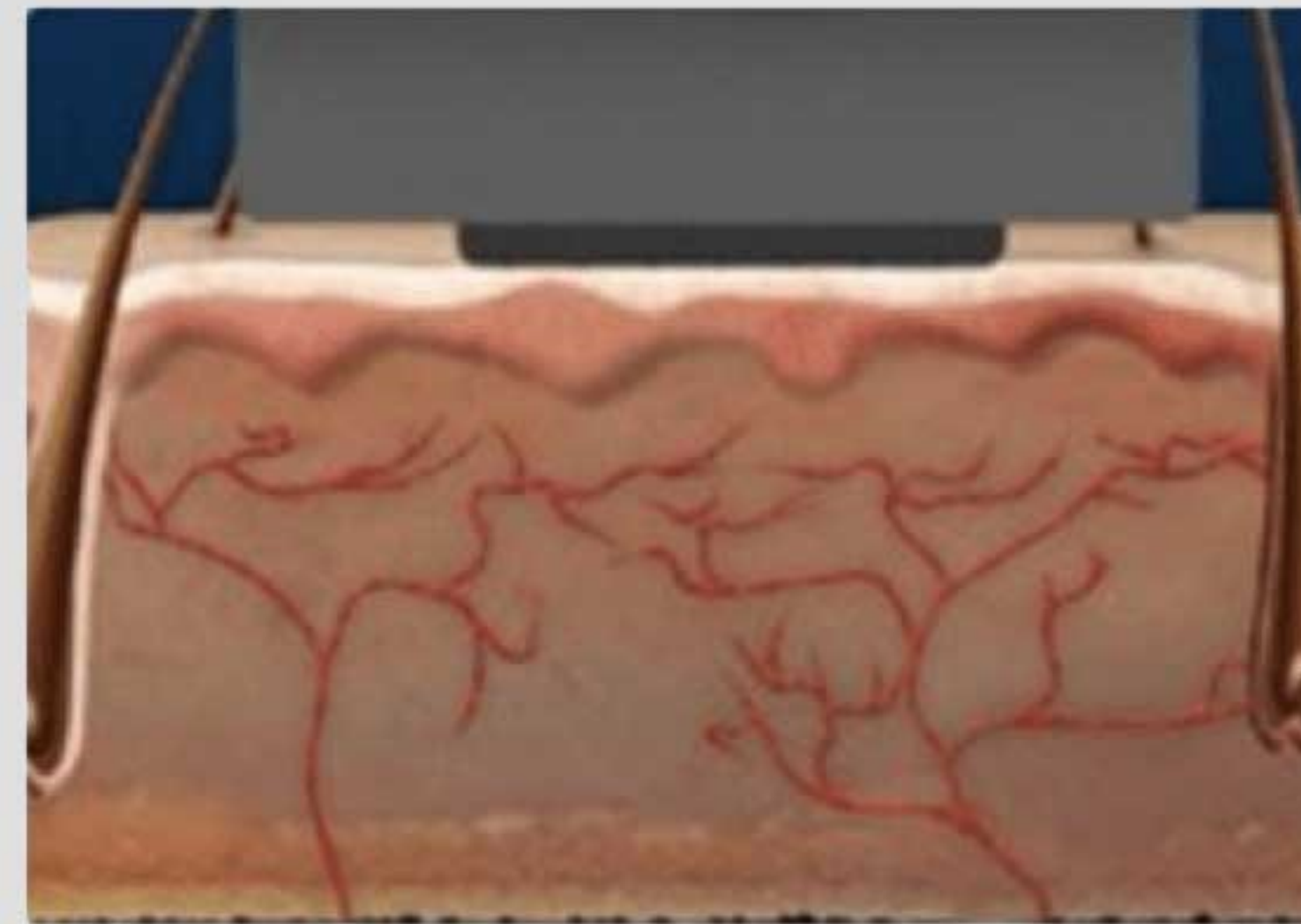


1step
Intracell makes contact with the skin

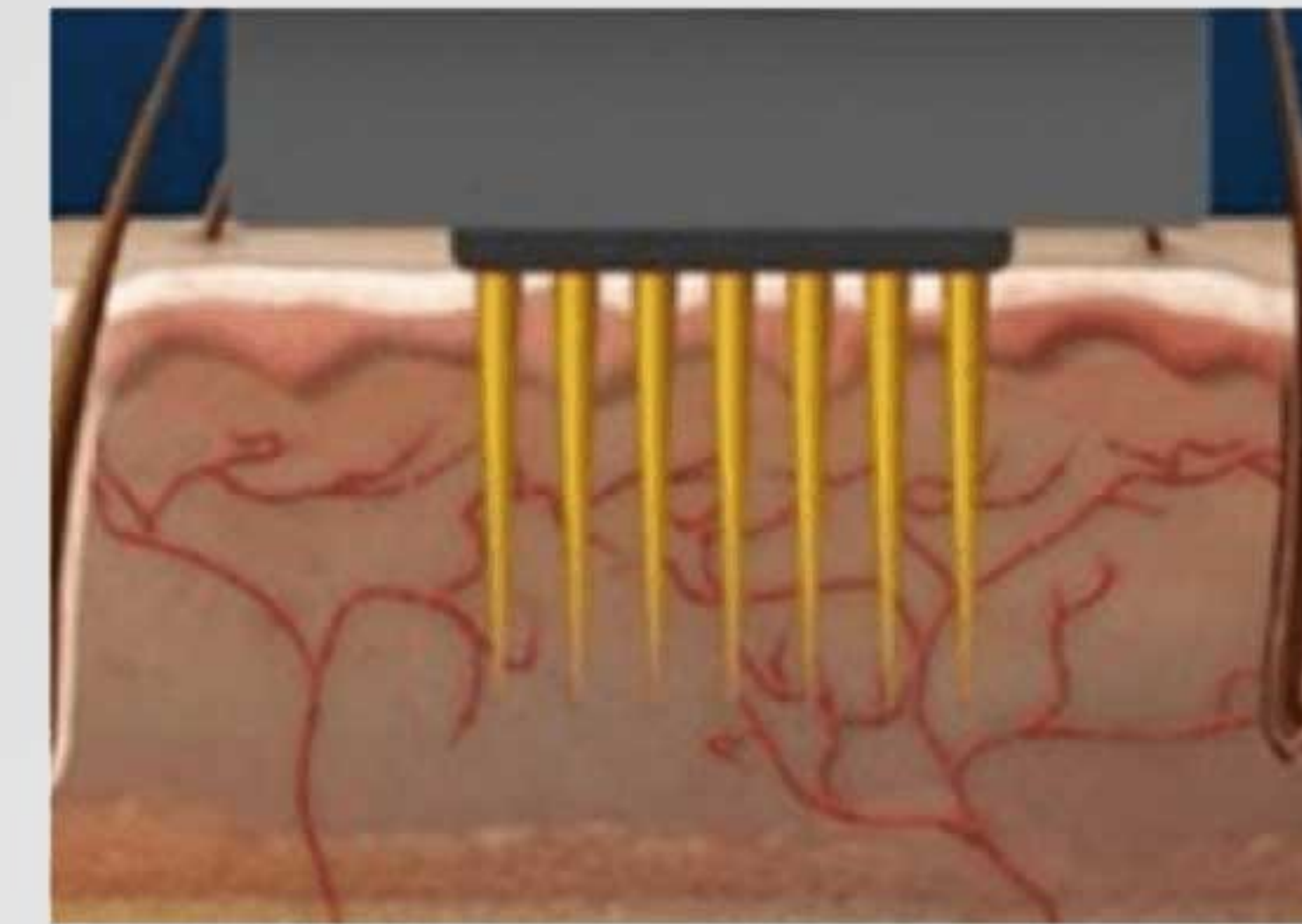


2step
Permeate to the desired depth

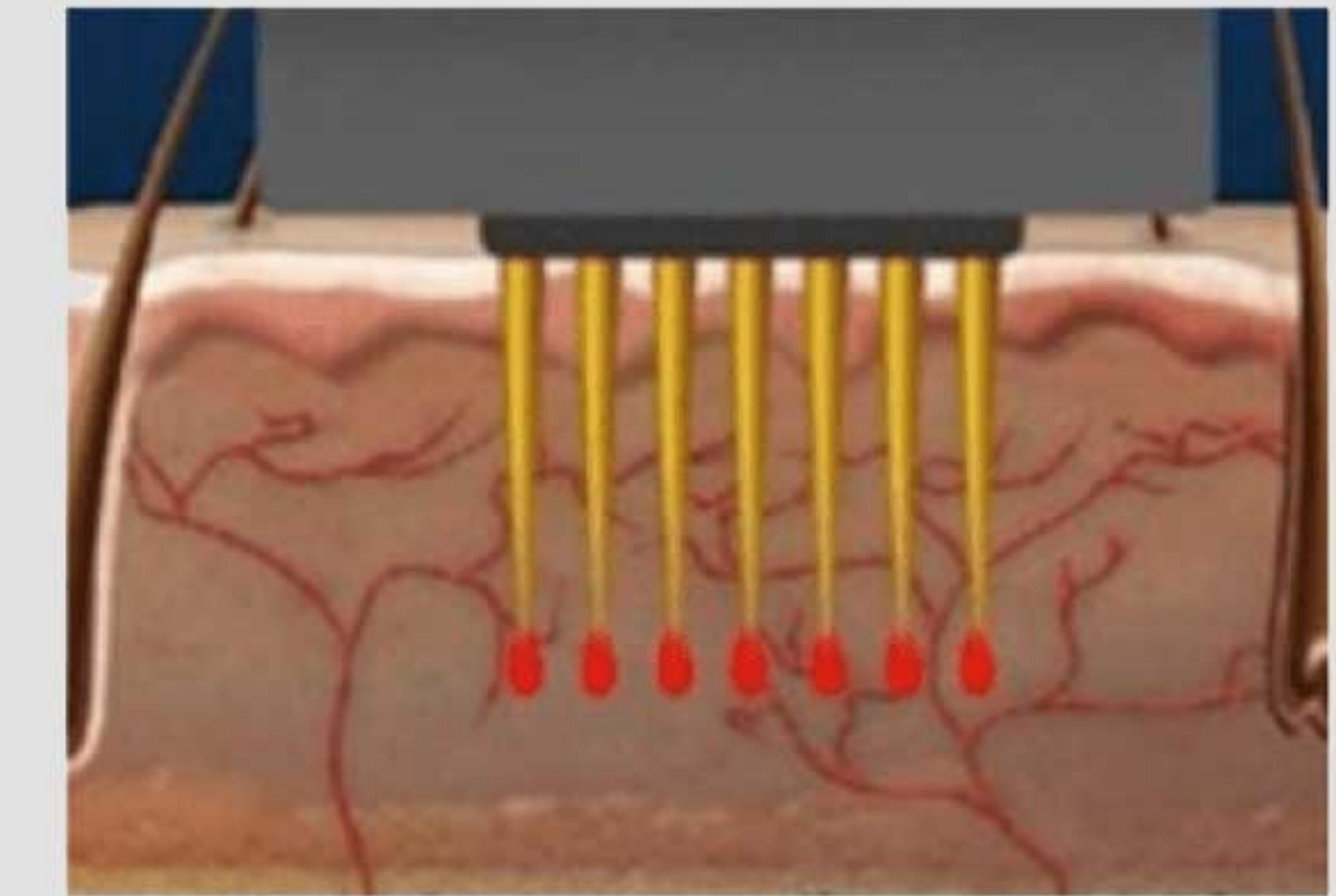


3step
Insonify through the gold needle

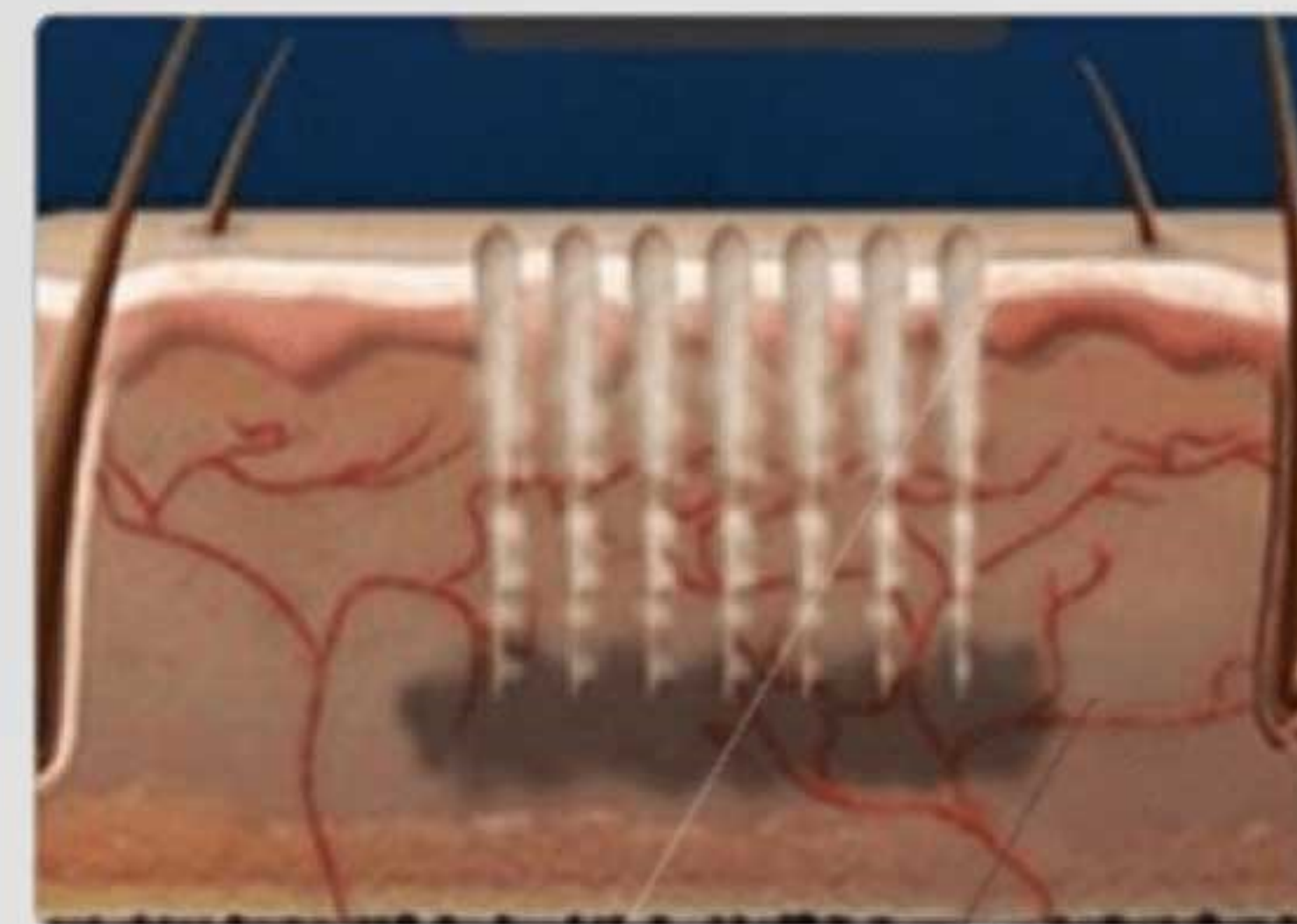

4step
Only the collagens at the point of the needle are destroyed

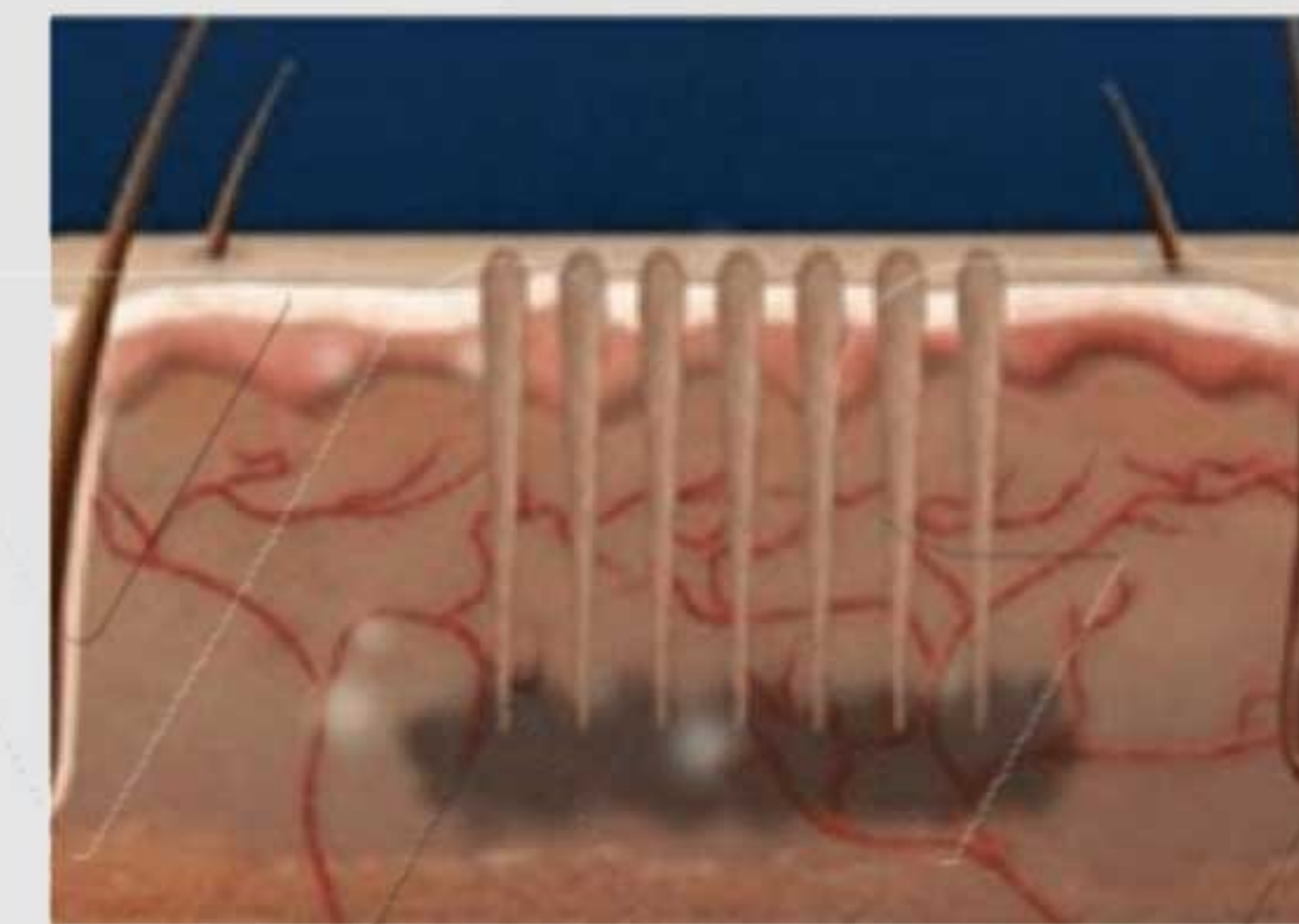


5step
Skin regeneration

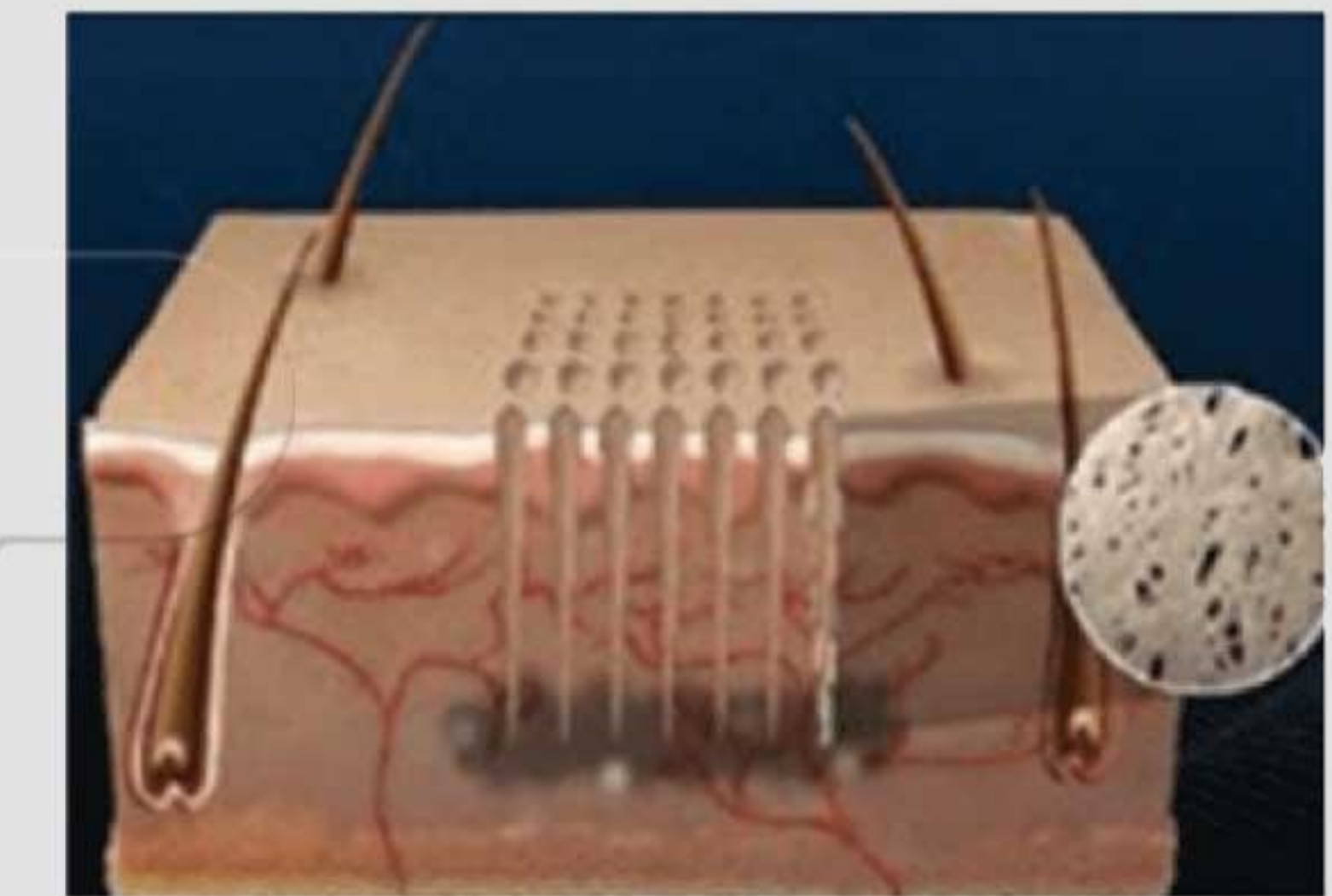

6step
Skin regeneration continues for the next 70 days


RF Tightening Lift (bipolar)

Radiofrequency skin tightening, or RF treatment, is a non-invasive procedure that some doctors do in their office. The doctor will use a wand that is connected to a machine that creates radio waves. They will pass the wand over your body to apply radio waves to your skin.

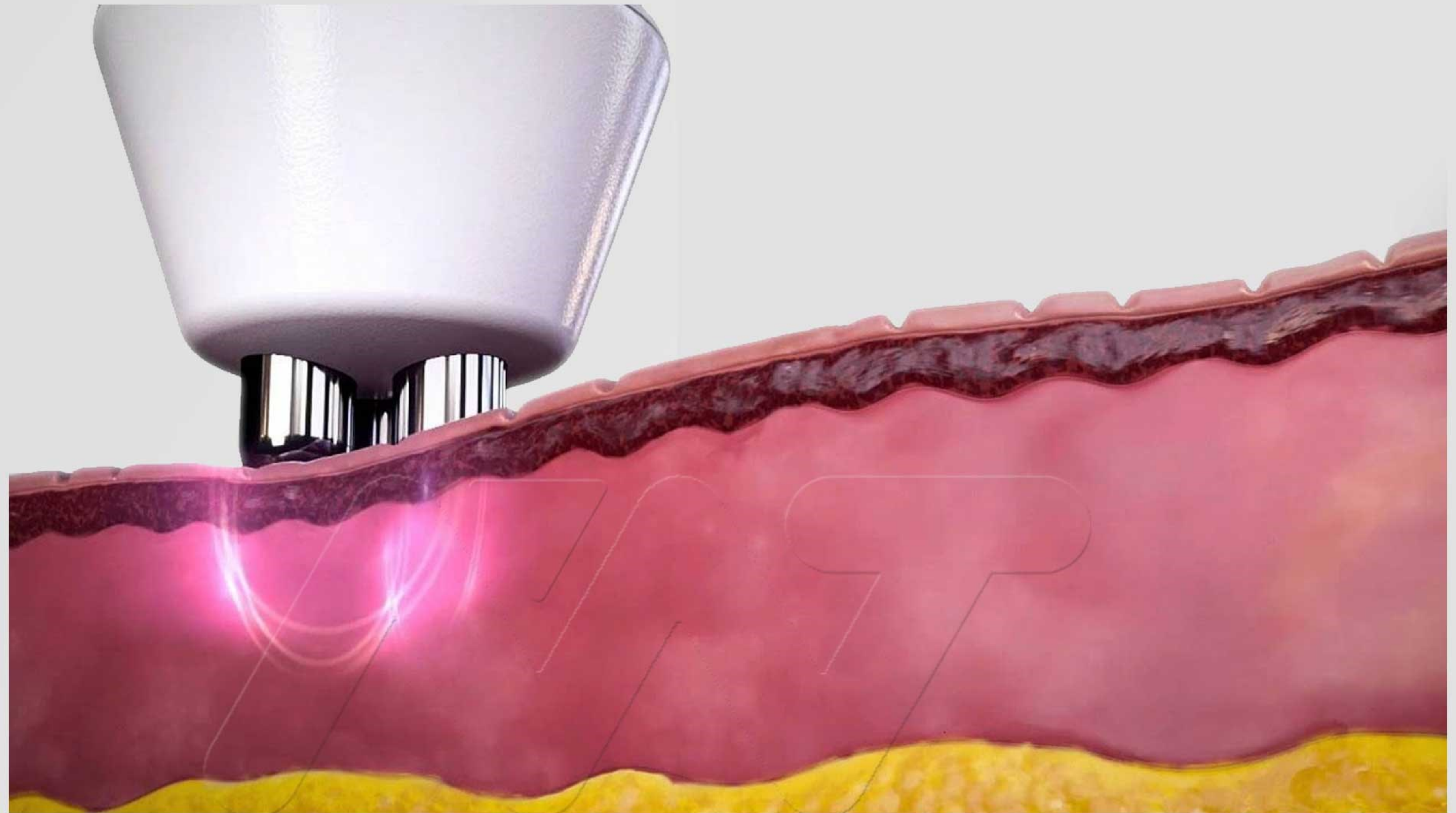

Ultrasound Face & Eye

An Ultrasound Facial creates high speed oscillations and uses low-frequency sound waves to leave the skin smooth, fresh, and perfectly prepared for product application. The treatment is stimulating, soothing, moisturizing, non-invasive and a safe method of treatment for all skin types.

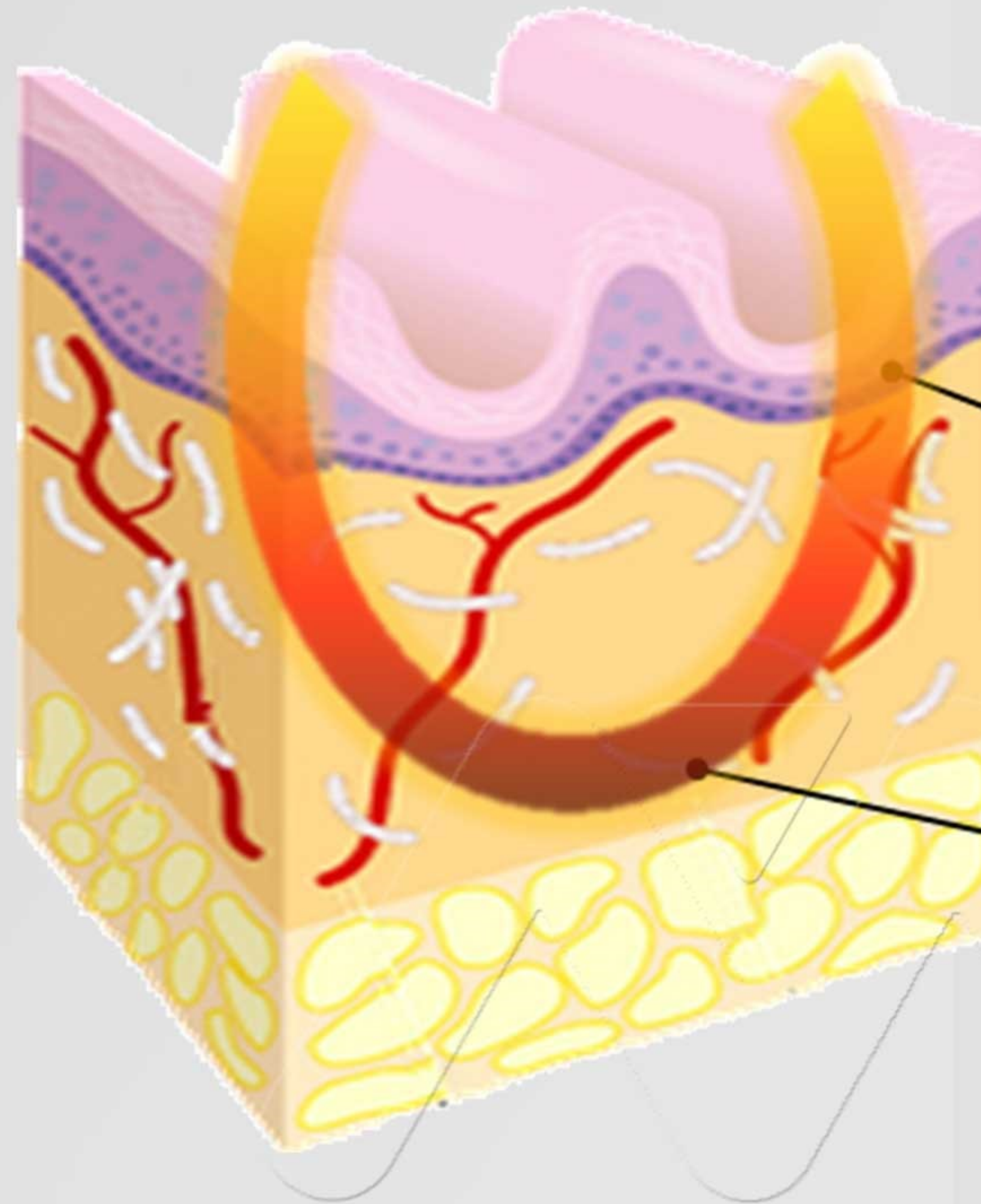


RF Energy penetrates into dermal tissues up to depth of 8mm

Dermal tissues are heated to approx. 42 degrees

Iontophoresis

Iontophoresis is a process of transdermal drug delivery by use of a voltage gradient on the skin. Molecules are transported across the stratum corneum by electrophoresis and electroosmosis and the electric field can also increase the permeability of the skin.

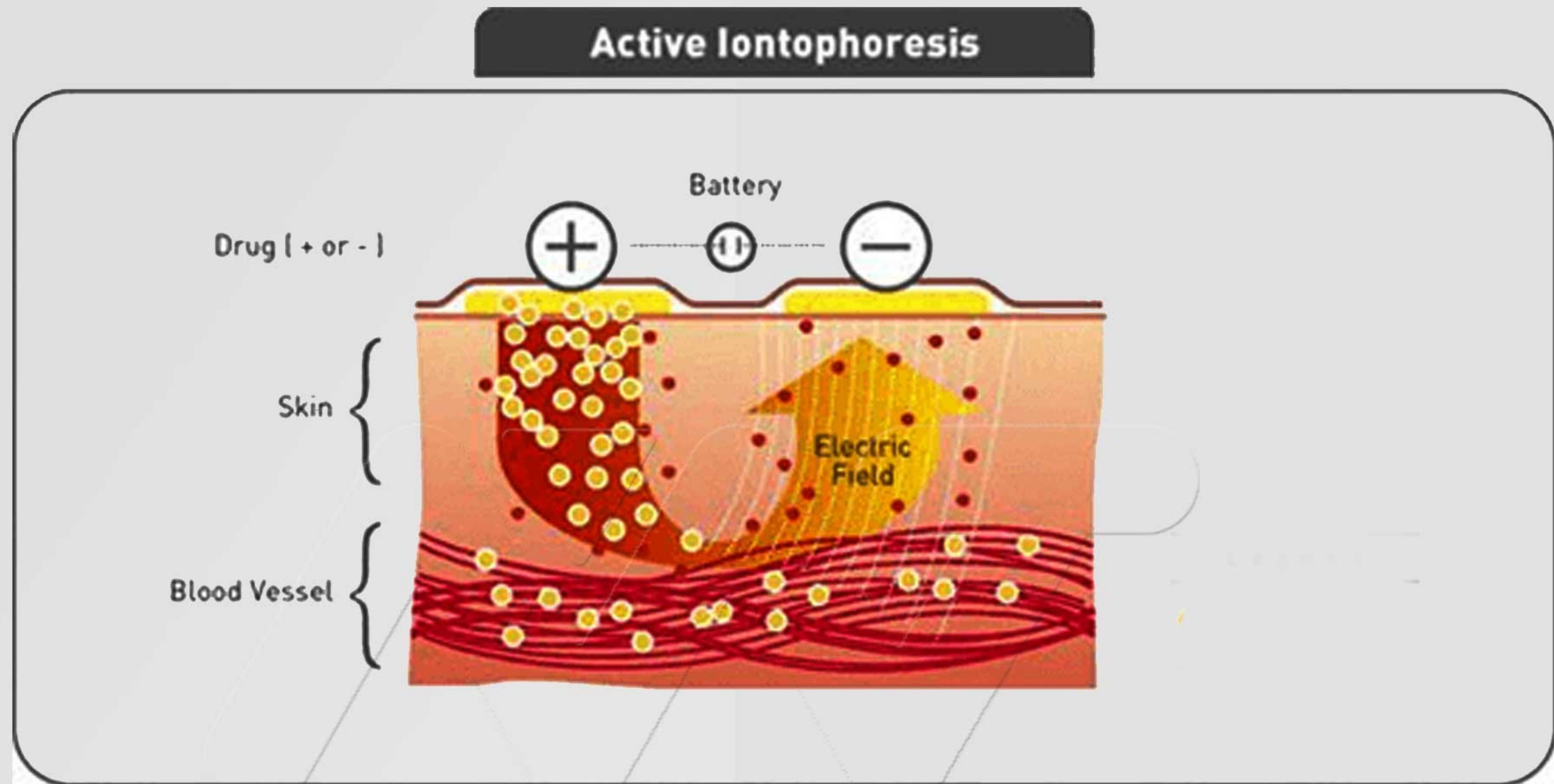


Cavitation

Cavitation is a non-invasive fat reduction treatment that uses ultrasound technology to reduce fat cells in targeted parts of the body. It is the preferred option for anyone who does not want to undergo extreme options such as liposuction, as it does not involve any needles or surgery.




How does ultrasonic cavitation work?

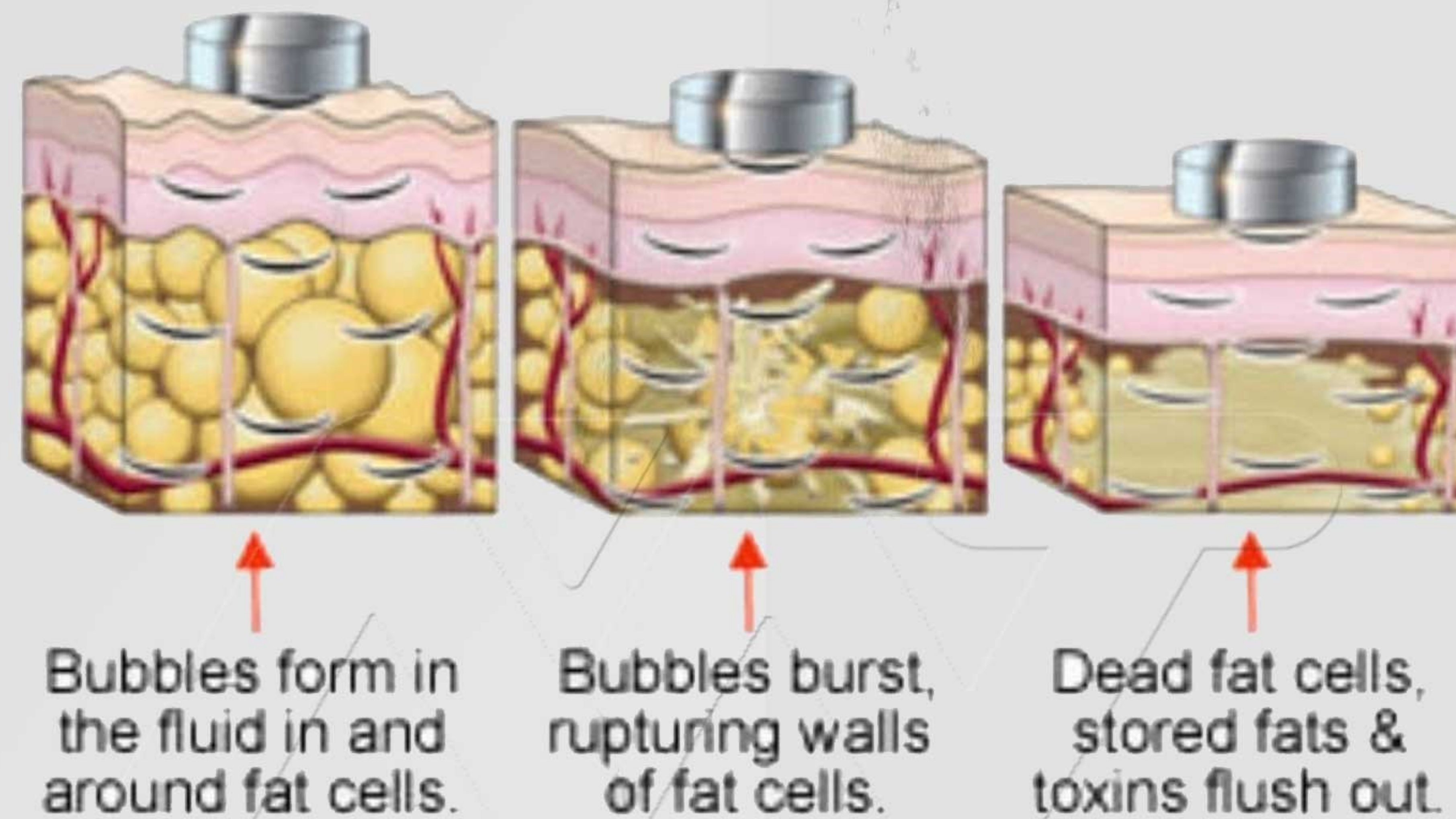

RF Lipolysis


Radiofrequency lipolysis is a procedure which uses a radio frequency to heat the fat without an instrument making physical contact with the patient.

The applicator heats from a distance of one centimeter from the skin.

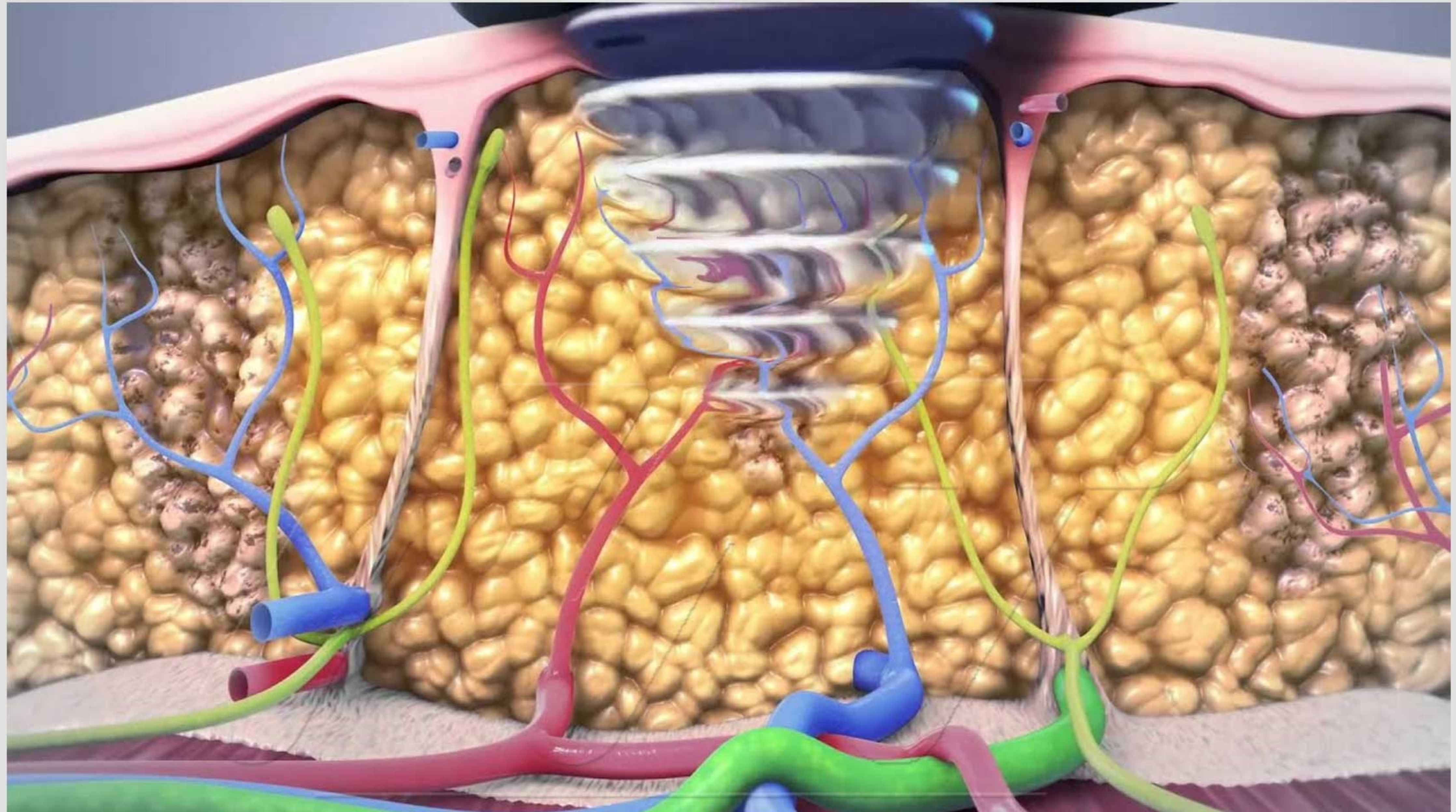

Roller Vacuum

Vacuum pressure and roller speed help in skin tightening, cellulite removal, and wrinkle removal. The use of different technologies plus different rollers and hand pieces cover the whole body treatments.

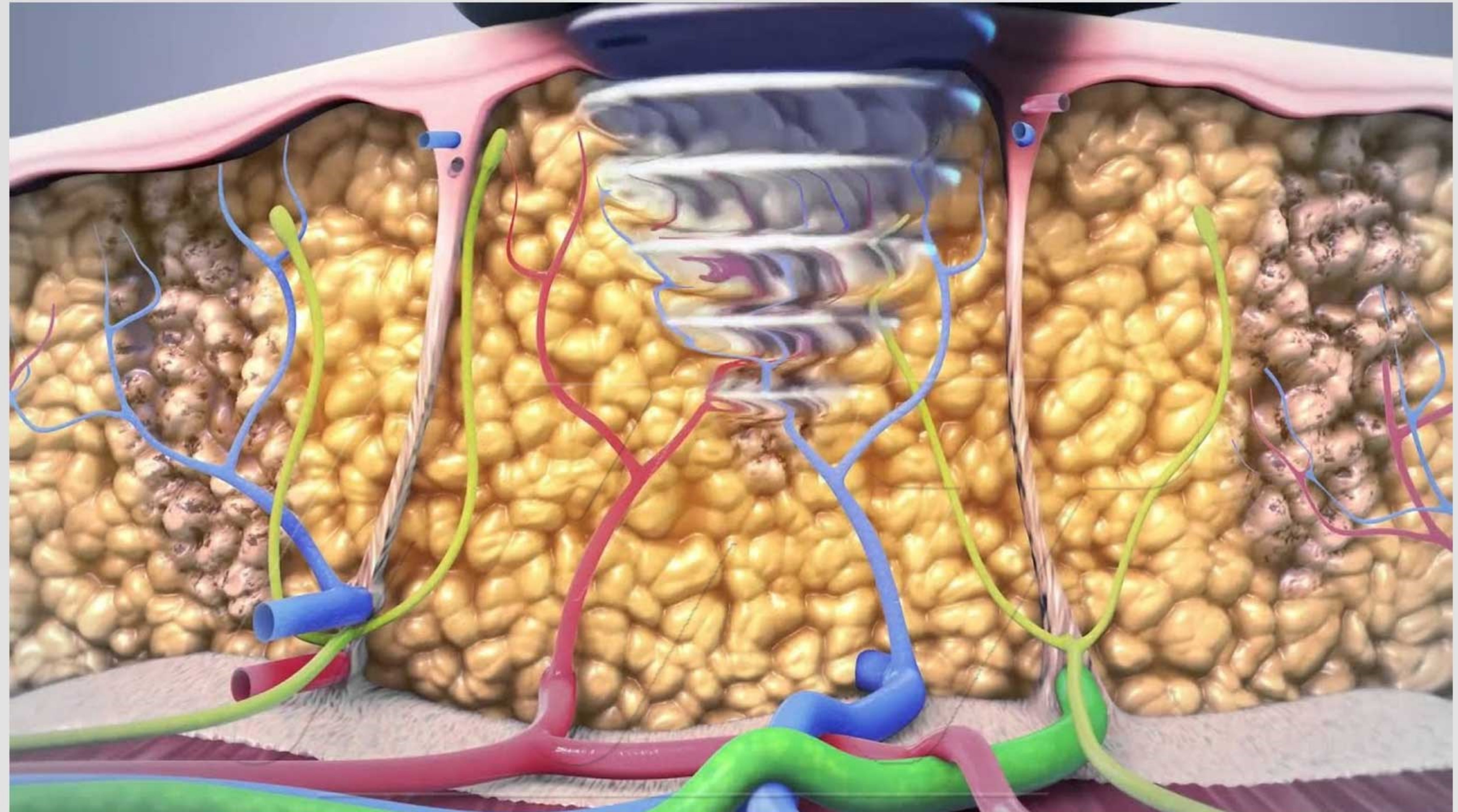



MULTITECH

021-22321127

09120783350

Mashhad, km 12, Asian Road - Science and Technology Park

www.faragard.co

fanavari.tech@gmail.com